


CLEANING AND DISINFECTION MACHINE TOPCLEAN 60S

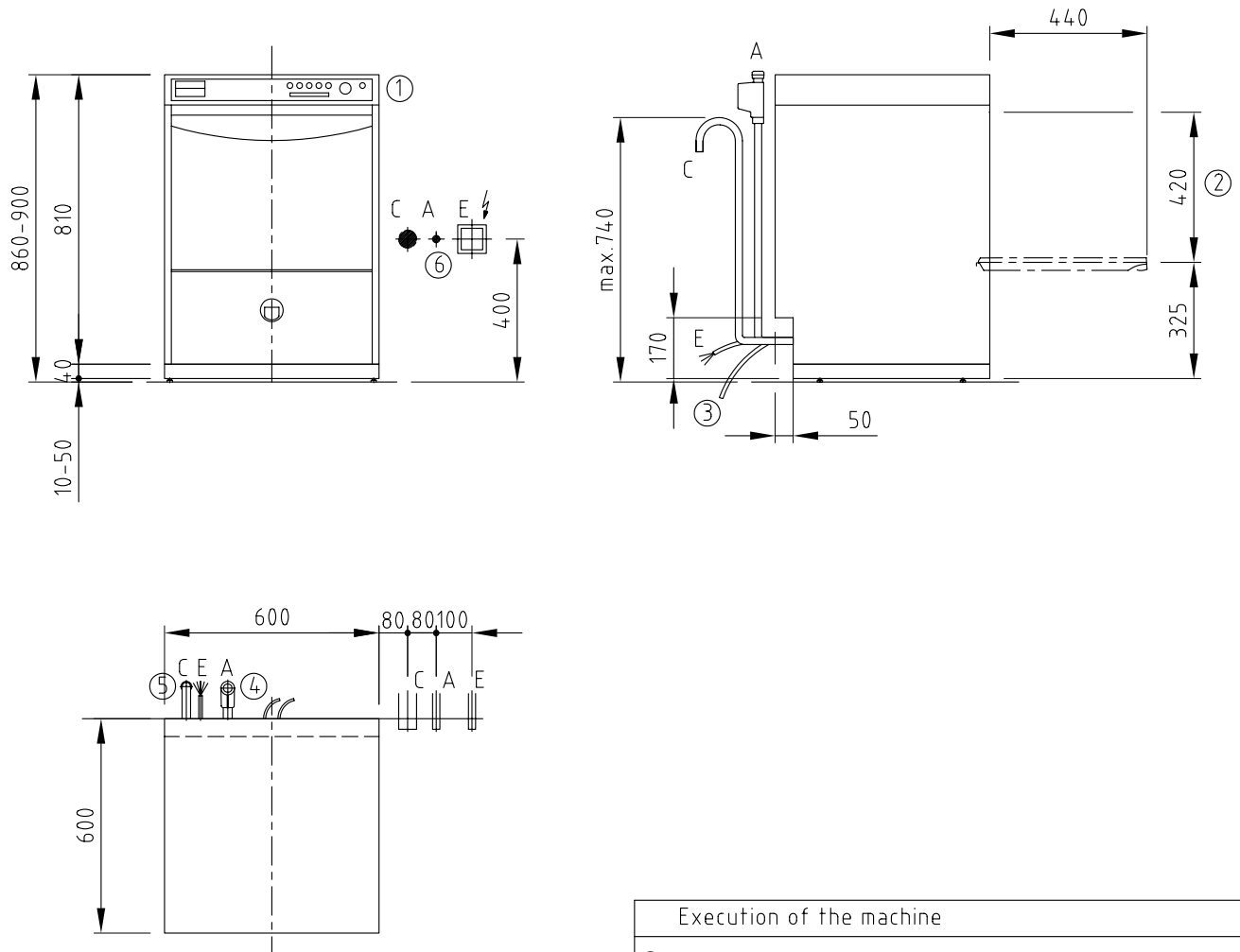

Execution of the machine	
①	Sensor touch panel
②	Entrance height
③	Rinse aid / Detergent supply line (option) with suction device
④	Connection kits possible at left or right side
⑤	Waste water hose inside the automat laid as syphon. Length of connection kit from exit of automat 1,4 m
⑥	Recommended local connecting points (possible at right or left hand side)
Subject to technical changes!	

A Sanitary installation			E El-connection		
Nominal value	Throughput	Water hardness	Feeding temperature	Total connections	Nominal current
Thread R 3/4"	5 l/min	max. 3°GH	10 ° C	12,19 kW	19,52 Amp.
C Drain			45-60 ° C	9,19 kW	15,19 Amp.
Drain pump with hose 3/4"			Voltage 3NPE 400V~50Hz		
Locally available:	Freshwater feeding DN20. Regulating tap and fine filter $\leq 25 \mu\text{m}$ Minimum Flow pressure 0,6 in front of solenoid valve Maximum pressure 5 bar Drain DN50 with slope		Locally available: Protection has to be executed according to total connections value and VDE 0100. Automat to be included in equipotential bonding. Main switch and main fuses locally		
Ind.	Révision	Date	Name	Order-No: Drawing-No: S00-P03164 4 Reference:	
			Scale: 1:20 Drawn: 29.04.08 Name:ej Checked: 29.04.08 Name:ej		